



Full XPE™ Cup

Acetabular Hip System



Full XPE™ Cup


Ease of Handling

- A proprietary locking mechanism allows all cup options to be inserted with the same cup positioner



PMMA Spacer

- Helps to centralize the cup and create a uniform 2.5 mm cement thickness intended to evenly distribute the load transfer to the cement and bone

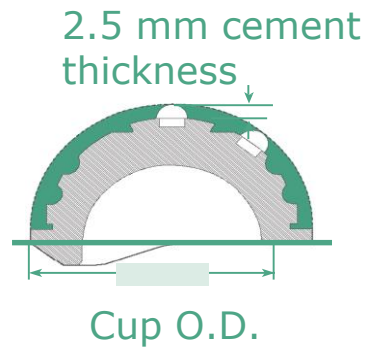

Horizontal & Vertical Groove

- Designed to enhance torque resistance and cement-implant fixation

Pressurizing Flange

- Designed to allow for optimal cement penetration and to help prevent cement extrusion.


Full XPE Cup and Femoral Head Compatibility Chart



Cup Size (mm)	Head Size (mm)
42	26
44, 46	28
48, 50	32
52~62	36

X-ray Marker

- Metal ring for X-ray image identification

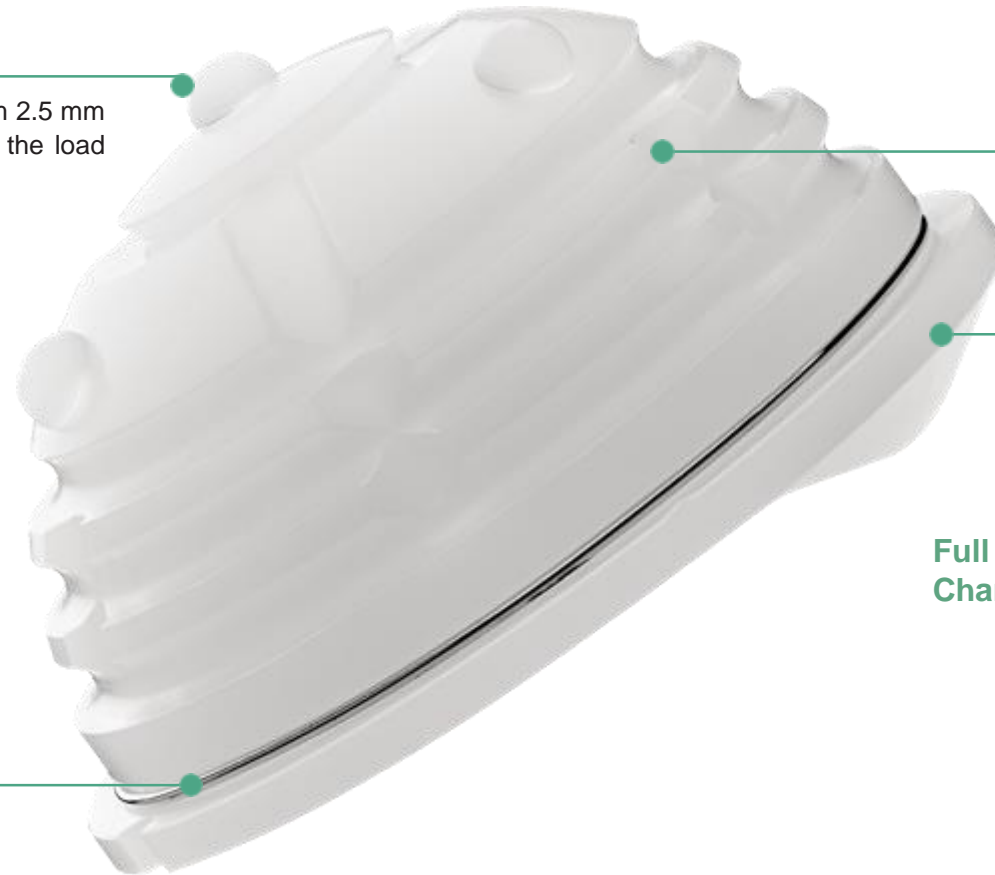

THANK YOU